

Extraction of THC, THCA and Carboxy-THC from Oral Fluid Using ISOLUTE® SLE+ after Collection with the Oral-Eze® Collection Device Prior to GC/MS Analysis

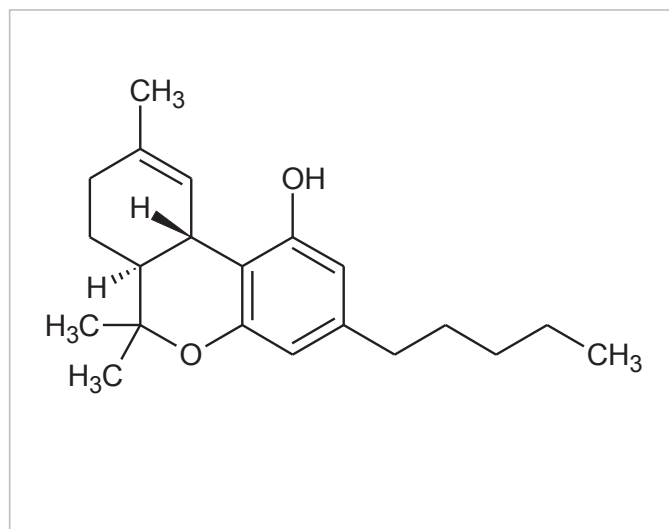


Figure 1. Structure of Δ^9 -THC (tetrahydrocannabinol)

Introduction

This application note describes the extraction of THC, THCA and Carboxy-THC from oral fluid matrix collected using the Oral-Eze® collection devices, prior to GC/MS analysis.

ISOLUTE® SLE+ Supported Liquid Extraction plates and columns offer an efficient alternative to traditional liquid-liquid extraction (LLE) for bio-analytical sample preparation, providing high analyte recoveries, no emulsion formation, and significantly reduced sample preparation.

This application note describes an effective and efficient ISOLUTE SLE+ protocol optimized for the 400 μ L sample capacity column format. The simple sample preparation procedure delivers clean extracts and excellent analyte recoveries.

Analytes

THC, THCA, THC-COOH and THC-d₃, THC-COOH-d₃ as internal standards

Sample Preparation Procedure

- Sample pre-treatment:** Following collection, add 4% ammonium hydroxide (aq) (10 μ L) to each collection device (see notes for buffer preparation instructions).
- Format:** **ISOLUTE® SLE+ 400 μ L Sample Volume Columns, part number 820-0055-B**
- Sample loading:** Load 300 μ L of the pre-treated oral fluid onto the column and apply a pulse of vacuum or positive pressure (3–5 seconds) to initiate flow. Allow the sample to absorb for 5 minutes.
- Analyte Extraction:** Apply dichloromethane/isopropanol, (95/5, v/v, 1 mL) and allow to flow under gravity for 5 minutes. Apply a further aliquot of DCM/IPA, (95/5, v/v, 1 mL) and allow to flow for another 5 minutes under gravity. Apply vacuum or positive pressure (5–10 seconds) to complete elution.
- Post Elution & Reconstitution:** Dry the extract in a stream of air or nitrogen using a SPE Dry (40 °C, 20 to 40 L/min) or TurboVap (1.0 bar at 40 °C for 40 mins).
- Upon dryness, reconstitute with 50 μ L ethyl acetate and 25 μ L MTBSTFA:TBDMCS 99:1 and vortex for 20 seconds. Transfer to a high recovery glass vial. Place in a heating block set to 70 °C, for 25 minutes. Remove vial from the block and allow cooling.

GC Conditions

Instrument:	Agilent 7890A with QuickSwap
Column:	Phenomenex Zebron ZB-Semivolatiles, 30 m x 0.25 mm ID x 0.25 µm
Carrier	Helium 1.2 mL/min (constant flow)
Inlet:	250 °C, Splitless, purge flow: 50 mL/min at 1.0 min
Injection:	2 µL
Wash solvents:	Ethyl acetate
Oven:	Initial temperature 100 °C Ramp 100 °C/min to 280 °C, hold for 10.5 minutes Ramp 100 °C/min to 330 °C, hold for 0.5 minutes
Post run:	Backflush for 2.4 minutes (3 void volumes)
Transfer Line:	280 °C

MS Conditions

Instrument:	Agilent 5975C
Source:	230 °C
Quadrupole:	150 °C
MSD mode:	SIM

SIM Parameters

Table 1. Ions acquired in the Selected Ion Monitoring (SIM) mode

SIM Group	Analyte	Target (Quant) Ion	1 st Qual Ion	2 nd Qual Ion
1	THC-d3	374	431	348
1	THC	371	428	345
2	THCA	530	631	455
3	THC-COOH-d3	416	518	575
3	THC-COOH	413	515	572

Results

The optimized ISOLUTE® SLE+ protocol demonstrated analyte recoveries ranging from 103–109% as shown in **Figure 2**. RSDs were below 10% for all analytes.

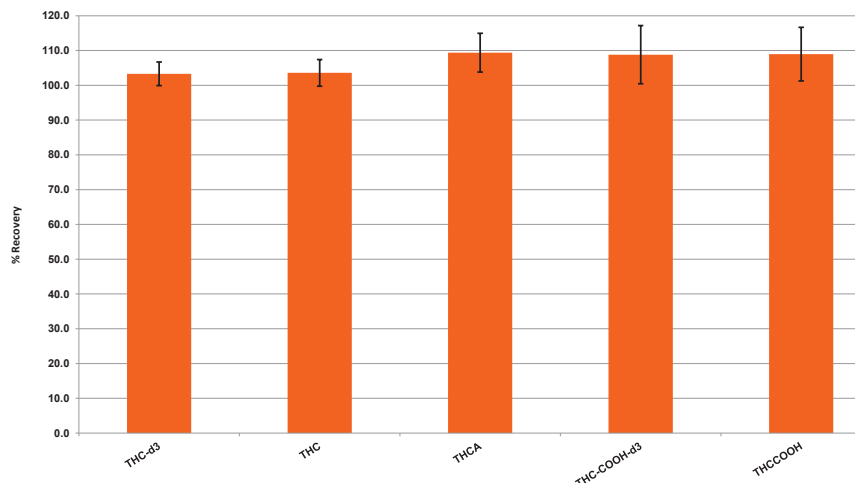


Figure 2. Typical extraction % recoveries (n=7) using the ISOLUTE® SLE+ protocol, at a concentration of 1 µg/mL of oral fluid.

Additional Information

If a non-chlorinated elution solvent is preferred, MTBE (methyl-tert-butyl-ether) is a suitable substitute solvent.

4% ammonium hydroxide is prepared from concentrated stock (28–30%) by adding 200 µL in 5 mL HPLC grade water.

Pre-treating to basic conditions accommodates multiple basic drugs of abuse in a single assay alongside THC, THCA and THC-COOH, if required.

Ordering Information

Part Number	Description	Quantity
820-0055-B	ISOLUTE® SLE+ 400 µL Supported Liquid Extraction Columns	50
PPM-48	Biotage® PRESSURE+ 48 Positive Pressure Manifold	1
SD-9600-DHS-EU	Biotage® SPE Dry Sample Concentrator System 220/240 V	1
SD-9600-DHS-NA	Biotage® SPE Dry Sample Concentrator System 100/120 V	1
C103198	TurboVap® LV, 100/120V	1
C103199	TurboVap® LV, 220/240V	1

For the latest application notes visit www.biotage.com

EUROPE

Main Office: +46 18 565900
 Toll Free: +800 18 565710
 Fax: +46 18 591922
 Order Tel: +46 18 565710
 Order Fax: +46 18 565705
order@biotage.com
 Support Tel: +46 18 56 59 11
 Support Fax: +46 18 56 57 11
eu-1-pointsupport@biotage.com

NORTH & LATIN AMERICA

Main Office: +1 704 654 4900
 Toll Free: +1 800 446 4752
 Fax: +1 704 654 4917
 Order Tel: +1 704 654 4900
 Order Fax: +1 434 296 8217
ordermailbox@biotage.com
 Support Tel: +1 800 446 4752
 Outside US: +1 704 654 4900
us-1-pointsupport@biotage.com

JAPAN

Tel: +81 3 5627 3123
 Fax: +81 3 5627 3121
jp_order@biotage.com
jp-1-pointsupport@biotage.com

CHINA

Tel: +86 21 2898 6655
 Fax: +86 21 2898 6153
cn_order@biotage.com
cn-1-pointsupport@biotage.com

To locate a distributor,
 please visit our website at
www.biotage.com

Part Number: AN825

© 2014 Biotage. All rights reserved. No material may be reproduced or published without the written permission of Biotage. Information in this document is subject to change without notice and does not represent any commitment from Biotage. E&OE. Product and company names mentioned herein may be trademarks or registered trademarks and/or service marks of their respective owners, and are used only for explanation and to the owners' benefit, without intent to infringe.

For more information visit www.biotage.com.